


Cordis[®]

OPTA[®] PRO

0.035" PTA Dilatation Catheter

Excellent crossability and
versatility to treat routine,
or challenging cases.

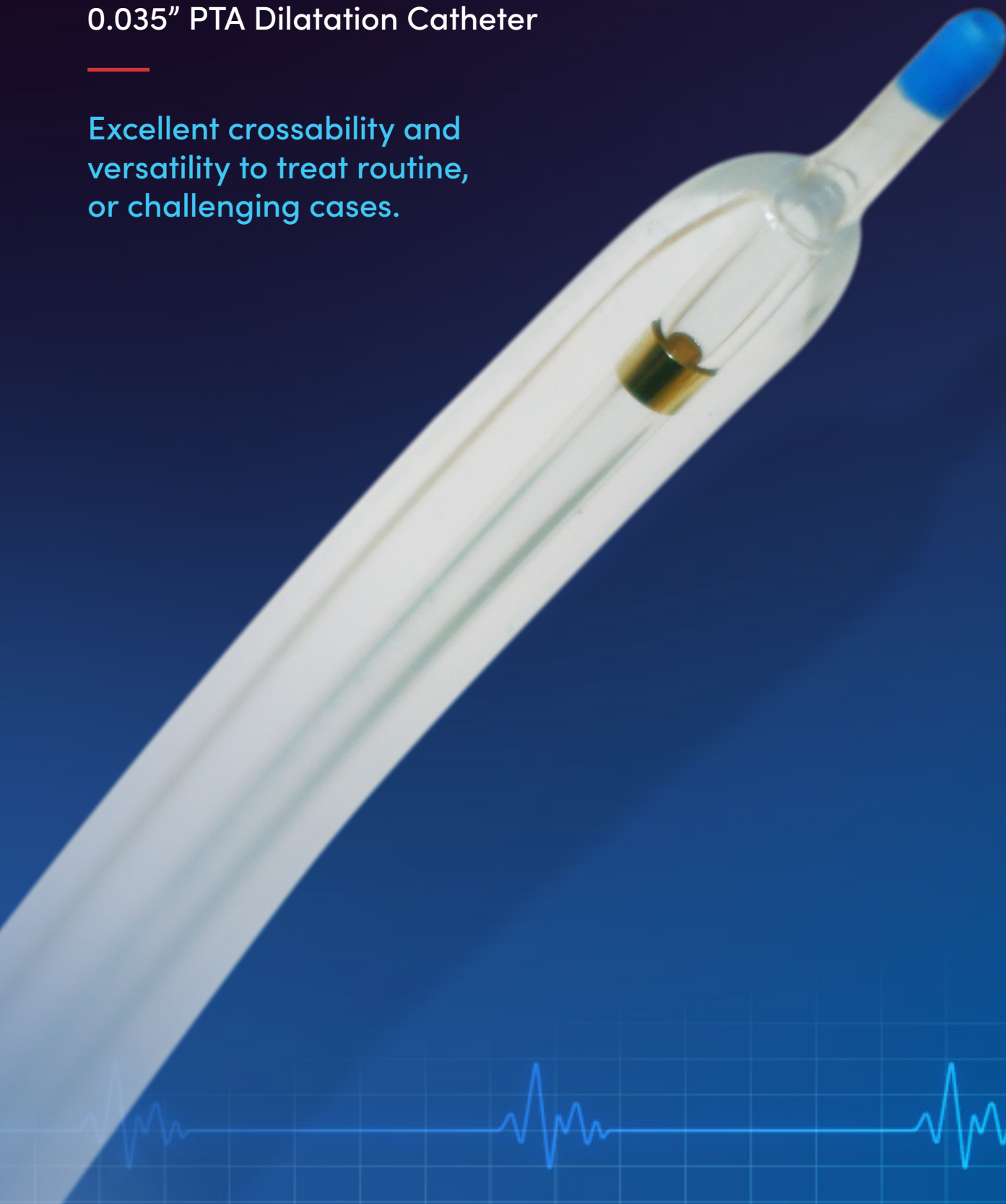

OPTA[®] PRO

0.035" PTA Dilatation Catheter

ADVANCED CROSSABILITY

- Flexible, tapered tip to ease entry into peripheral arteries to facilitate the crossing of tight stenoses.
- MDX Silicone coated Distal (30cm) shaft and guidewire lumen to facilitate smooth delivery and wire handling.
- DURALYN™ combines flexibility and strength for improved crossability and dilation of long lesions.
- Puncture and abrasion resistant balloon material.

VERSATILITY


- Rated Burst Pressures
 - » 10 atm on 3mm to 10mm sizes with nominal at 6 atms.
 - » 6 atm on 12mm with nominal at 4 atm.
- DURALYN™ High Performance Balloon Material provides controlled compliance for predictable balloon sizing.
- Two gold, radiopaque markers facilitate placement of the OPTA[®] PRO PTA Catheter for effective dilatation.

ACCESS

- Low Profile:
 - » Wide array of 5F compatible sizes.
 - » Excellent trackability and pushability to facilitate lesion access.

EFFECTIVE

- Effective in treatment of stenoses in peripheral arteries.
- Effective in the treatment of obstructive lesions of native or synthetic arteriovenous dialysis fistulae.

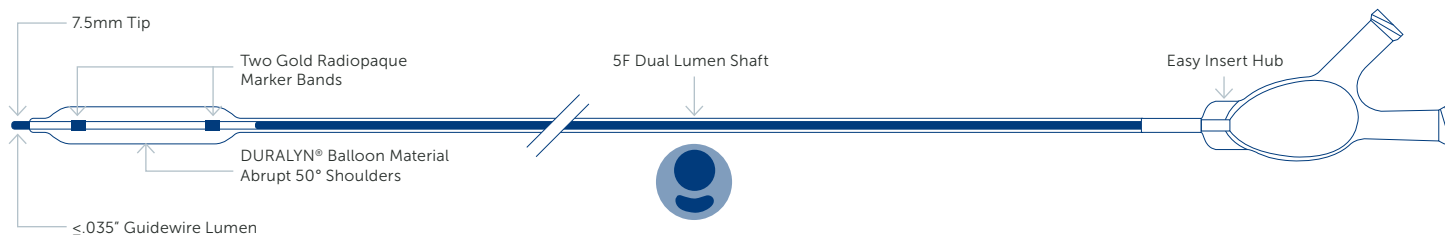

FEATURES

- 5F
- Dual lumen shaft
- Flexible Tapered Tip
- Optimal rated burst pressure: 10 atm (3mm - 10mm) with a nominal pressure of 6 atm
- Extended product range
- One piece hub
- Inner and outer (30cm distal) MDX coating

OPTA[®] PRO

0.035" PTA Dilatation Catheter

ORDERING INFORMATION

Packaging: 1 unit/Box

Balloon Diameter (mm)	Balloon Length (mm)	Rated Burst Pressure (atm)	Sheath Fit (F)	Shaft Length (80cm)	Shaft Length (110cm)	Shaft Length (135cm)
3	40	10	5	4193040S	4193040L	
3	80	10	5		4193080L	4193080X
4	20	10	5	4194020S		4194020X
4	40	10	5	4194040S		
4	80	10	5	4194080S	4194080L	
4	100	10	5	4194000S	4194000L	
5	20	10	5	4195020S		
5	30	10	5	4195030S		
5	40	10	5	4195040S		4195040X
5	100	10	5	4195000S	4195000L	
6	20	10	5	4196020S	4196020L	
6	40	10	5	4196040S	4196040L	4196040X
6	60	10	5		4196060L	
6	80	10	6		4196080L	
6	100	10	6		4196000L	
7	20	10	5	4197020S	4197020L	
7	30	10	5	4197030S		
7	40	10	5	4197040S	4197040L	
7	60	10	6	4197060S	4197060L	
7	80	10	6		4197080L	
7	100	10	6	4197000S		
8	20	10	6	4198020S	4198020L	
8	30	10	6	4198030S		
8	40	10	6	4198040S	4198040L	
8	60	10	6	4198060S		
8	80	10	6	4198080S	4198080L	
8	100	10	7	4198000S		
9	20	10	6	4199020S	4199020L	
9	40	10	6	4199040S		
10	20	10	6	4190020S	4190020L	
10	30	10	6	4190030S		
10	40	10	6	4190040S	4190040L	4190040X
10	60	10	7	4190060S		
12	20	6	7	4192020S	4192020L	
12	30	6	7	4192030S		
12	40	6	7	4192040S	4192040L	

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician.

Important information: Prior to use, refer to the instructions for use supplied with this device for indications, contraindications, side effects, suggested procedure, warnings and precautions.

CORDIS, Cordis LOGO, OPTA, and DURALYN are trademarks of Cordis and may be registered in the US and/or in other countries.

© 2021 Cordis. All Rights Reserved. 100590400 09/2021